



TIPS:

Your suggestions and comments for Costway are really important to us!
We sincerely solicit you to go back to our shop and leave a good rating in just a simple click. It would be quite encouraging if you could kindly do so like below:

★★★★★ Great products so far


February 24, 2018

Great product so far. Fast delivery, easy setup, and working without any issues.

With your inspiring rating, Costway will be more consistent to offer you
EASY SHOPPING EXPERIENCE, GOOD PRODUCTS and EFFICIENT SERVICE!

US office:Fontana

UK office:Ipswich

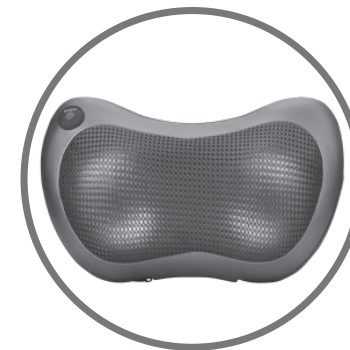


USER'S MANUAL

**Massage Pillow
EP23474US**



In order to keep the massager in the best condition at any time, it is recommended not to use more than 30 minutes for every time, and let the machine rest for 30 minutes after each use.

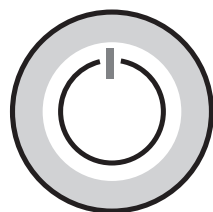

Switch diagram

Note: 1. the working time of product should not exceed 30 minutes.
2. The machine is equipped with overheating protection, the motor will be automatically shut down when overheating. The machine can work again after it is cooled

Using Method

Pressing the first time: the automatic mode the heating function will be started, after 1 minute of forward moving, it automatically switches to reverse moving. After 1 minute of the reversal moving, it automatically switches to forward moving, and then the cycle is repeated again and again.

Pressing again: after 1 minutes of reverse moving, it then automatically switches to forward moving. after 1 minutes of forward moving, it automatically switches to reverse moving, and the cycle is repeated again and again.

Pressing again to turn off the heating function.

And then pressing again: massage pillow stops working.

◆ **The using method for the speed switch (limited to the machine with double button and speed control functions)**

The default speed adjustment is the high-speed mode, press the speed switch once for the medium-speed mode, and then press once again for the low-speed mode, and further press once to return to the high-speed mode, the cycle is repeated again and again.



Recommended to drink 500cc warm water before and after using the product to help detoxification, and promote metabolism.

Product attention and maintenance

In order to keep the massage pillow to be in the best condition, please note that the following questions and only use the product in accordance with the functional description of the manual.

1. Do not massage the surface of the body with trauma or ulceration.
2. To prevent this product from overloading, do not stand or pressurize on it.
3. When you start using the product, please do not use it for a long time, and you should gradually increase the time duration of use after you get used to the product.
4. The time of use should not be more than 30 minutes each time.

5. During the course of the use, when the product 's temperature is too high or you smell the bad gas from it, please stop using it and turn off the power.
6. The pillow can produce surface heating, the user who is heat-insensitive must be careful when using it.
7. Use a soft cotton cloth and a mild detergent to clean it.
8. Do not use (or store) this product in wet or hot places.
9. This product can not be stamped, heavily thrown or pressurized.
10. Do not use this product if you are receiving medical treatment and any of the following conditions.
 - People with cancer
 - The body is equipped with a pacemaker
 - People with chronic or major illness
 - Pregnant women and children
 - Suffering from severe heart disease
 - Mental illness
 - The body is fitted with other metal utensils, such as prostheses
11. Do not energetically pull the power cord of the massage pillow.
12. Do not use this massage pillow in water or in a wet place. If it accidentally falls into the water, do not take the massage pillow immediately out of the water, you should first turn off the power and remove its power plug from the outlet.
13. Do not use the massage pillow near electric heaters and other heating sources.
14. Do not use accessories or parts supplied by the other companies on this massage pillow.
15. Do not open the massage pillow or replace the parts, or the warranty will be void, and the relevant after-sales service will be stopped.
16. For your personal safety and the life of this product, do not use this massage pillow with zipper opening. The Company will not be responsible for any consequences arising from the violation of maintenance and use rules.
17. Do not use this product within 1 hour before and after meals.
18. If the power cord or plug is damaged, please do not use further, it must be sent to the manufacturer for repairing or replacement to avoid danger.

Technical parameters

Product Name: massage pillow
 Rated power: 24W
 Output voltage: DC12V --- 2A
 Safety design:

Product model: EP23474
 Input voltage: AC100-240V 50/60Hz
 Rated time: 15 minutes